



URF TYPE **P5**






BREAKTHROUGH TECHNOLOGY EVERY ASPECT OF THE ENDOUR

From insertion to maintenance, Olympus incorporates its latest proprietary technology in the development of the URF TYPE P5 Flexible Ureteroscope where significant benefits are realized when performing uretero-renoscopy.

Access

Flexible angulation facilitates orientation inside the kidney.

- The 275° down angulation enables optimal visualization in the lower calyx.
- The 180° up angulation (with a small radius) is most appropriate for accessing the upper/middle calyces.
- The scope is uniquely designed to maintain excellent tip deflection even when small caliber laser fiber is inserted.

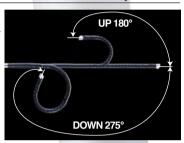

Visualization

Image quality has been significantly improved.

By incorporating a moiré-reduction filter into the URF TYPE P5 scope, Olympus has succeeded in further enhancing its industry superior opto-digital imaging performance. The moiré-reduction effect provides the user with an extremely clear picture.







Control

One hand operation, improved torque control, detachable LG

- While angulating the scope with one hand, you can simultaneously operate a forceps or open/close the irrigation port with the other hand.
- A special new material has been designed into the outer sheath making tip rotation more responsive when twisting the control section.
- Detachable light guide (LG) capability allows the user to quickly switch back and forth the light guide between the cystoscope/semi-rigid ureteroscope and the URF TYPE P5.



ADVANCES OLOGIC PROCESS

Maintenance

Remarkable durability & leakage testing

- A special new material used in the outer sheath minimizes crushing.
- Thicker rubber in the bending section resists potential puncturing.
- Leakage testing allows punctures to be detected early and helps to ensure a longer service life for the scope.










while maintaining optimal image quality and treatment

capabilities. Because of the slim 5.3 Fr. tip diameter,

insertability is significantly enhanced.

Evolutiontip

n Accessories

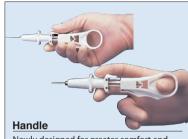


Newly designed for greater comfort and efficiency. Can be used with two styles of grip.



FG-51D/52D/55D (Length:1150mm) Disposable Grasping Forceps (Basket Type)

Features a basket with large openings. Designed to not be crushed in narrow canals.

*FG-51D/55D: 3 Fr., 4 wires; FG-52D: 2.4 Fr., 3 wires

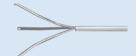

FG-54D (3Fr., Length:1150mm) Disposable Grasping Forceps (Three Nail Type)

Specially designed for the difficult task of removing small calculi from the pelvis



FB-56D-1(3Fr., Length:1150mm) Biopsy Forceps

Autoclavable forceps that can be inserted into a 3 Fr. (1.2 mm) channel.

A-----

MD-82 Guide Wire

0.035 inch diameter and flexibility on both sides provides smoother insertion.



ST-U1 Urinary Duct Sheath Dilator

Its slim design (11.7 Fr.) lets it fit easily on the URFTYPE P5. Also, up to 0.035 inch guide wire can be used for smoother insertion



WA05991A Instrument Tray



n Specifications

Optical System	Field of View	90°
	Direction of View	0° (Forward Viewing)
	Depth of Field	2 ~ 50mm
Insertion Tube	Distal End Outer Diameter	5.3 Fr. (bullet shape)
	Insertion Tube Outer Diameter	8.4 Fr. (ø2.8mm)
	Working Length	700mm
Instrument Channel	Inner Channel Diameter	3.6 Fr. (ø1.2mm)
Bending Section	Angulation Range	Up 180°/Down 275°
Total Length		1050mm

n Standard Set


OES Uretero-Reno Fiberscope URF TYPE P5 — 1	MH-507 Channel-Opening Cleaning Brush ————
MAJ-579 Biopsy Valve — 10	BW-15SH Channel Cleaning Brush —————
MAJ-891 Forceps/Irrigation Plug 1	BW-7B Channel Cleaning Brush —————
MB-156 FTO Cap 1	

n Endourology Products



EndoEYE Ureteroscope

With an advanced CCD built into the distal tip, this innovative video scope will astonish you with its superb, true-to-life image quality. The light guide cable and remote control switches are integrated within the body of the scope.

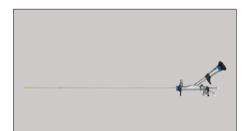

OES Pro Ureteroscope

The compact 50,000-image-fiber bundle provides excellent visualization. Larger, brighter, more consistently high-quality images enable detailed observation and effective treatment. Expansive irrigation capability increases efficiency. Uniquely designed insertion portion and revised tip profile provides atraumatic insertion. Autoclavable.



OES Pro Percutaneous Nephroscope

This new percutaneous nephroscope offers unsurpassed image quality. Produces excellent flow capability, ensuring the clearest field of view possible and allowing you to concentrate on your procedures. A variety of forceps are available. Ergonomically designed. Autoclavable.



CYF TYPE V2/VA2 VISERA CYSTO-NEPHRO VIDEOSCOPE

The flexible video cystoscope incorporates a high-resolution, multicolor CCD chip at the distal end of the scope which produces a large, bright, and consistent image. The Evolutiontip design enhances insertability without compromising image quality and treatment capability. A suction function is also available with the CYF-VA2.



CYF-5/5A OES CYSTO NEPHRO FIBERSCOPE

Compatible with electrosurgical procedures.

One-handed control of the aspiration function makes it easy to keep the field of view clear. The large ø2.4mm (7.2 Fr.) channel allows perfusion with 6 Fr. instruments. Evolutiontip design.



LUS-2 Ultrasonic Lithotriptor

Highly effective stone fragmentation and suction using ultrasonic power. Integrated generator and pump (two-in-one system) improves both procedural and cost efficiency.

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.



OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan OLYMPUS WINTER & IBE GMBH Kuehstrasse 61, 22045 Hamburg, Germany OLYMPUS MEDICAL SYSTEMS EUROPA GMBH Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany OLYMPUS SURGICAL & INDUSTRIAL AMERICA INC. One Corporate Drive, Orangeburg, New York 10982, U.S.A (100 Bible Lagoon Drive, Suite 390 Miami, Florida 33126-2087, U.S.A KEYMED LTD.

Kewmed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

OLYMPUS SINGAPORE PTE LTD.
4918 River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373
OLYMPUS HONG KONG AND CHINA LIMITED
Room 1520-1527, Soon Centre, 15 Gartish Road, Taimshatsul, Kowloon, Hong Kong
Room, 1202-161, Soon Centre, 15 Gartish Road, Taimshatsul, Kowloon, Hong Kong
Room, 1202-161, Soon Centre, 15 Gartish Road, Taimshatsul, Kowloon, Hong Kong
LTD.
Room, 1202-161, Tower, At 2, languamena, Marenus Chaoyang Olstrict Beijing, 100222, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071 Moscow, Melaya Kaluzirakaya 19, 154, 11, 12, Russia
OLYMPUS AUSTRALIA PTY LTD
31 Gilby Road, Mount Waverley, VIC, 3143 A ustralia